

User Guide: Managing Webhooks in Quartzzy

[Overview](#)

[What are Webhooks and their uses?](#)

[How do Webhooks work?](#)

[Access webhook configurations](#)

[Create a new webhook configuration](#)

[Modify an existing webhook configuration](#)

[Enable or Disable an existing webhook configuration](#)

Overview

What are Webhooks and their uses?

Webhooks can be used to build integrations between Quartzzy and another software tool used in your lab.

- These integrations are outbound in nature – which means that they help you share information from Quartzzy to this other software tool.
- These integrations will help you save cost and increase efficiency by automating tasks in a lab thereby cutting down on manual work and reducing human errors.

Examples of such outbound integration use cases include:

- Communication and Task management
 - Keep track of your Quartzzy order status updates in Slack/Microsoft Teams
 - Create a task in Asana or a reminder in Google calendar and prepare ahead for upcoming deliveries
- Reporting/Analytics
 - Automatically keep track of your monthly spend on Quartzzy orders in Google Sheets and save time otherwise spent copying-pasting information from Quartzzy
- Inventory Management
 - Sync information about Quartzzy inventory items to your favorite database tools such as Airtable, Notion, Google Sheets, etc.

How do Webhooks work?

Webhooks allow you to receive events as HTTP POST requests to a server when events occur in Quartzzy. These events can be used to design these outbound integration use cases described in the previous section.


See [this Wikipedia article](#) for more information on the basics of a webhook. To learn more about best practices around Quartzzy's Webhooks, please read [this introductory article](#) in our developer documentation.

The upcoming sections will go into more detail about how to:

- Access webhook configurations
- Create a new webhook configuration
- Modify an existing webhook configuration
- Enable or Disable an existing webhook configuration

Access webhook configurations

To access existing webhook configurations or create a new configuration, [navigate to your Lab Settings in Quartzly](#) and click on the tab named **Webhooks**. Here is an example below:

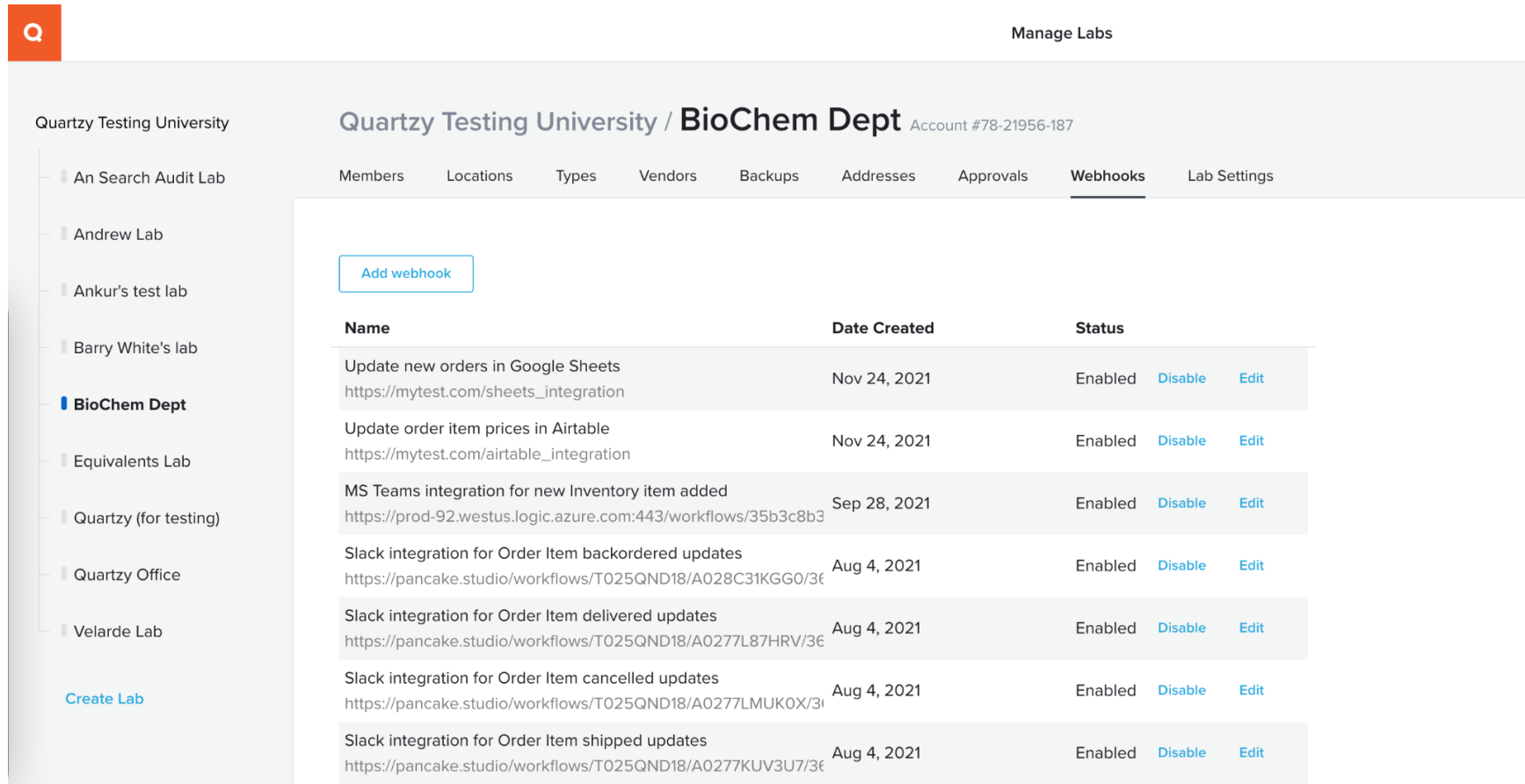


The screenshot displays the 'Manage Labs' interface in Quartzly. On the left, a sidebar lists various labs, with 'BioChem Dept' selected. The main content area shows the configuration for 'BioChem Dept' (Account #78-21956-187). A navigation bar at the top of the main area includes tabs for 'Members', 'Locations', 'Types', 'Vendors', 'Backups', 'Addresses', 'Approvals', 'Webhooks', and 'Lab Settings'. The 'Webhooks' tab is active, showing a table of existing webhooks. An 'Add webhook' button is located at the top left of the table area.

Name	Date Created	Status
Update new orders in Google Sheets https://mytest.com/sheets_integration	Nov 24, 2021	Enabled Disable Edit
Update order item prices in Airtable https://mytest.com/airtable_integration	Nov 24, 2021	Enabled Disable Edit
MS Teams integration for new Inventory item added https://prod-92.westus.logic.azure.com:443/workflows/35b3c8b3	Sep 28, 2021	Enabled Disable Edit
Slack integration for Order Item backordered updates https://pancake.studio/workflows/T025QND18/A028C31KGG0/36	Aug 4, 2021	Enabled Disable Edit
Slack integration for Order Item delivered updates https://pancake.studio/workflows/T025QND18/A0277L87HRV/36	Aug 4, 2021	Enabled Disable Edit
Slack integration for Order Item cancelled updates https://pancake.studio/workflows/T025QND18/A0277LMUK0X/36	Aug 4, 2021	Enabled Disable Edit
Slack integration for Order Item shipped updates https://pancake.studio/workflows/T025QND18/A0277KUV3U7/36	Aug 4, 2021	Enabled Disable Edit

Create a new webhook configuration

1. To create a new webhook configuration, click the **Add webhook** button.
2. In the new modal called **Add webhook**, enter the three mandatory fields
 - a. Name: The display name of your webhook configuration
 - b. URL: The URL that would receive the webhook
 - c. Event Types: The event(s) occurring in Quartzly for which you want to receive webhooks
3. Click **Save**
4. In the new modal called **New webhook**, you will notice a field displaying the Secret Key. Click the copy icon to copy this secret key in a safe location. Click **Close** after you are done copying the secret key.

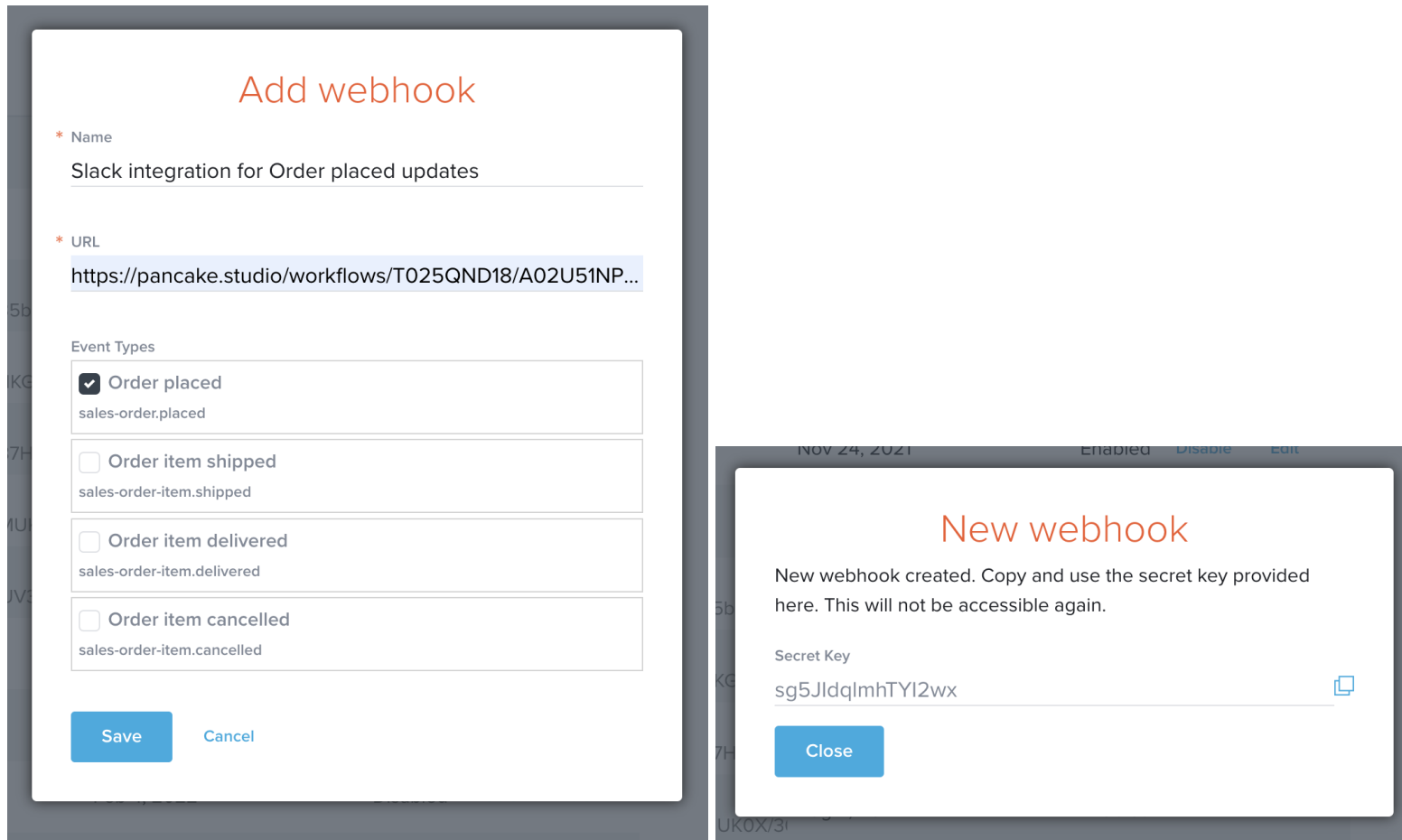

Modify an existing webhook configuration

1. To modify an existing webhook configuration, click the **Edit** link next to the webhook name in the table displaying all your webhook configurations.

2. In the new modal called **Edit webhook**, update the fields of interest.
3. Once you are done editing, click **Save**.

Edit webhook

* Name
Update order item prices in Airtable

* URL
https://mytest.com/airtable_integration

Event Types

- Order placed
sales-order.placed
- Order item shipped
sales-order-item.shipped
- Order item delivered
sales-order-item.delivered
- Order item cancelled
sales-order-item.cancelled

Save Cancel

Enable or Disable an existing webhook configuration

1. To enable or disable an existing webhook configuration, click the **Enable** or **Disable** link next to the webhook name in the table displaying all your webhook configurations.
2. Once you are done, you would notice the latest status of your webhook configuration in the **Status** column. You would also notice that on refreshing the page, the disabled webhook configurations appear at the bottom of the table as shown in the example below.



Quartz Testing University

Quartz Testing University / BioChem Dept Account #78-21956-187

- Members
- Locations
- Types
- Vendors
- Backups
- Addresses
- Approvals
- Webhooks**
- Lab Settings

[Add webhook](#)

Name	Date Created	Status
Slack integration for Order placed updates https://pancake.studio/workflows/T025QND18/A02U51NPR8T/36	Feb 23, 2022	Enabled Disable Edit
Update new orders in Google Sheets https://mytest.com/sheets_integration	Nov 24, 2021	Enabled Disable Edit
Update order item prices in Airtable https://mytest.com/airtable_integration	Nov 24, 2021	Enabled Disable Edit
MS Teams integration for new Inventory item added https://prod-92.westus.logic.azure.com:443/workflows/35b3c8b3	Sep 28, 2021	Enabled Disable Edit
Slack integration for Order Item backordered updates https://pancake.studio/workflows/T025QND18/A028C31KGG0/36	Aug 4, 2021	Enabled Disable Edit
Slack integration for Order Item delivered updates https://pancake.studio/workflows/T025QND18/A0277L87HRV/36	Aug 4, 2021	Enabled Disable Edit
Slack integration for Order Item cancelled updates https://pancake.studio/workflows/T025QND18/A0277LMUK0X/36	Aug 4, 2021	Enabled Disable Edit
Slack integration for Order Item shipped updates https://pancake.studio/workflows/T025QND18/A0277KUV3U/36	Aug 4, 2021	Enabled Disable Edit
Test 212121 google.com	Feb 8, 2022	Disabled Enable Edit
Webhook created at 2022-02-08 00:53:19 google.com	Feb 8, 2022	Disabled Enable Edit
New webhook 2 https://example.com/mywebhook	Feb 4, 2022	Disabled Enable Edit
Webhook created at 2022-02-01 21:37:16 https://example.com/quartz	Feb 1, 2022	Disabled Enable Edit

- An Search Audit Lab
- Andrew Lab
- Ankur's test lab
- Barry White's lab
- BioChem Dept**
- Equivalents Lab
- Quartz (for testing)
- Quartz Office
- Velarde Lab

[Create Lab](#)